





CUTTING EDGE
sense of innovation

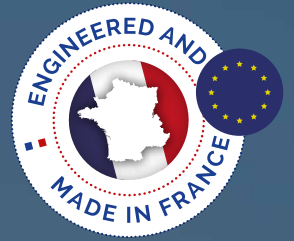

Synthesis[®]

Acrylic Monofocal

NEW

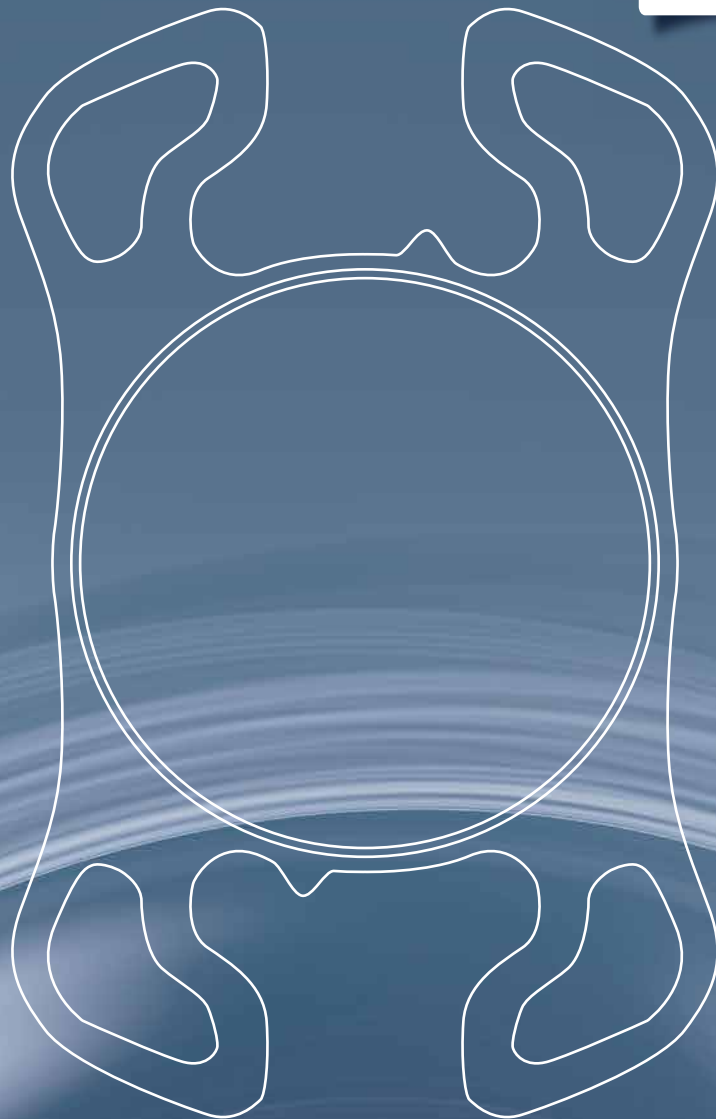

Synthesis[®]

Acrylic Monofocal

CUTTING-EDGE TECHNOLOGY IOL

Synthesis[®] is a new acrylic monofocal intraocular lens, developed and manufactured by CUTTING EDGE, expert in IOLs design.

Cutting Edge builds experienced and dedicated teams to design, prototype and manufacture high-quality Intraocular Lenses. Our state-of-the-art facility honours the quality and process sophistication of the medical device industry.

The proximity of our Research & Development team to our manufacturing operations is truly significant for building close links and in using a results-orientated approach.


Cutting Edge masters all manufacturing stages in house, from machining to optical and mechanical tests, through to packaging and international logistics.

Synthesis was engineered by a team who fully believe that the human eye is both fragile and precious and requires the utmost in precision, “savoir-faire” and faultless industrial quality.



Synthesis[®]

Acrylic Monofocal

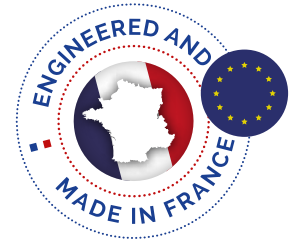

THE LATEST ALL-IN-ONE INNOVATIONS PLATFORM



SAFETY

- **Material**
 - Clear or blue-light-filtering option
- **Design**
 - 4 point fixation haptics
 - 360° C-Edge+ square edges
 - Orientation notch on optic edge
- **Truly MICS enhanced characteristics**

REPRODUCIBILITY

- **Optic**
 - 6 mm diameter optic
 - Anatomically friendly with 3 overall diameter
 - Aberration neutral Aspheric lens
 - Posterior shift features
- **Injection**
 - Latest injection system options

From -10.0 D to +40.0 D

- 1 D increment between -10.0 D to +10.0 D & +31.0 D to +40.0 D
- 0.5 D increment between +10.0 D to 30.0 D

OPTICAL QUALITY

- MTF curves exceed ISO 13485 norm
- UV and blue-light-filtering option



CUTTING EDGE

SAFETY

MATERIAL

- Acrylic copolymer successfully **used for more than 17 years** in ophthalmology
- **Ultraviolet absorption**
- **Blue-light-filtering** option (yellow)

DESIGN

- **4 point** fixation haptics
- **3 overall diameters** to fit the capsular bag:
 - 11 mm from -10.0 D to +15.0 D
 - 10.7 mm from +15.5 D to 22.0 D
 - 10.5 mm from +22.5 D to +40.0 D
- **Enhanced C-Edge+** $\leq 5\mu\text{m}$ 360° square edge

Synthesis[®]

Acrylic Monofocal



Overall diameter

- 11 mm
- 10,7 mm
- 10,5 mm

C-EDGE+ SCANNING ELECTRON MICROSCOPY IMAGES


MAGNIFICATION 50X



MAGNIFICATION 200X



MAGNIFICATION 500X


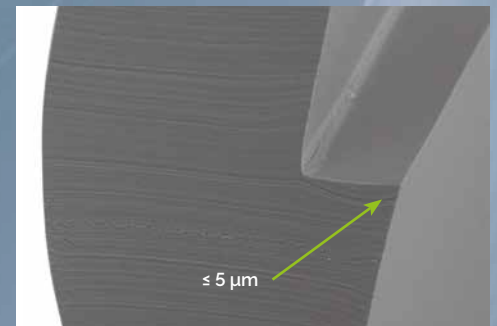


C-EDGE+

REPRODUCIBILITY

OPTIC


- **6 mm optic** across the power range
- **Orientation notch** at the optic edge
- **Posterior shift:** vaulted design for excellent capsular bag adhesion
- **Aberration neutral**
Aspheric profile
- **Homogeneous surface** across the power range for consistent injection forces regardless of IOL diopter.

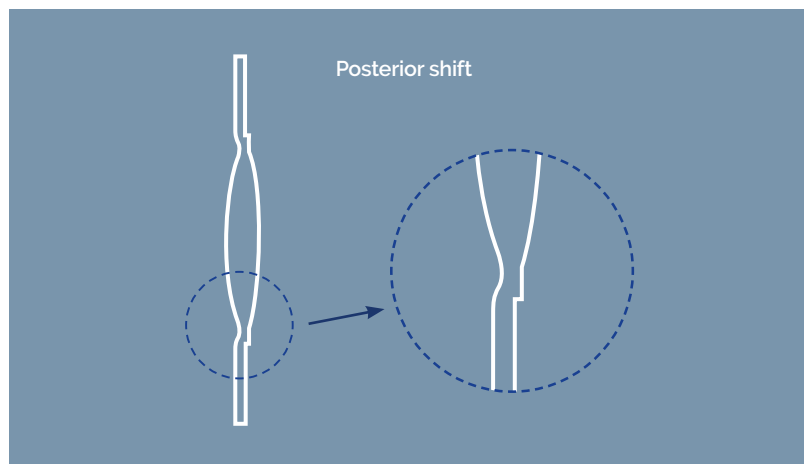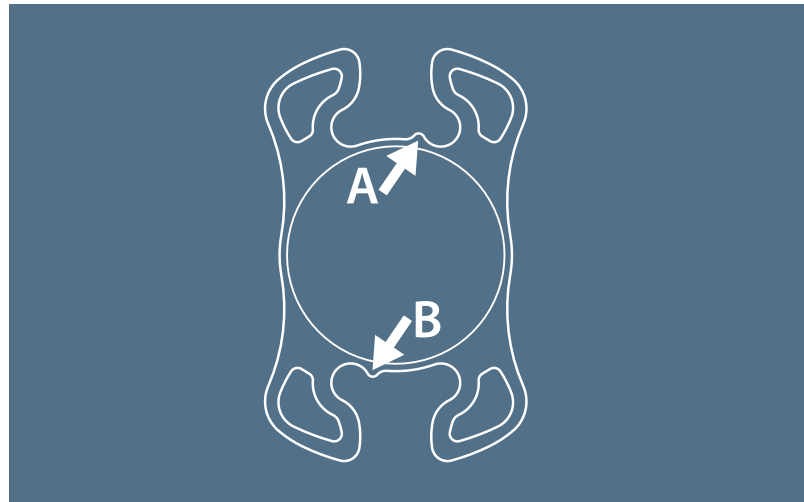

Synthesis[®]
Acrylic Monofocal



CE
0197

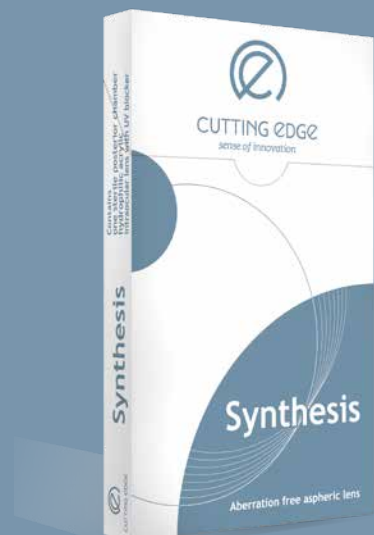


OPTICAL QUALITY


MTF (MODULATION TRANSFER FUNCTION)

Synthesis exceeds ISO 11979-2 standard requirement for optical properties and test methods (0.43 at 100 cycle/mm) at all light conditions and pupil diameter.

**SYNTHESIS® MTF
IN DAY LIGHT CONDITIONS
(3 mm pupil diameter)**

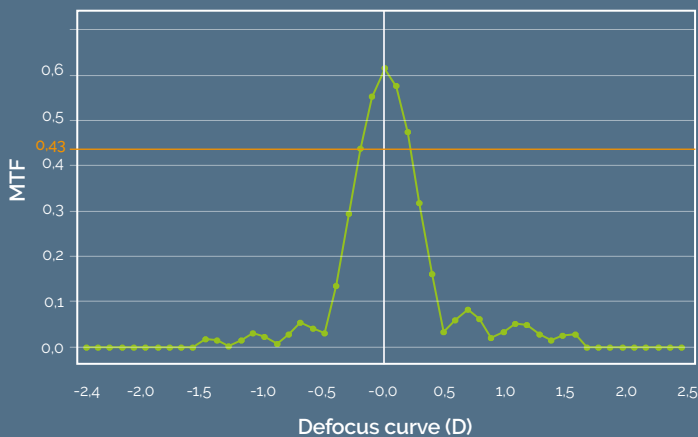


**SYNTHESIS® MTF
IN DIM LIGHT CONDITIONS
(4.5 mm pupil diameter)**

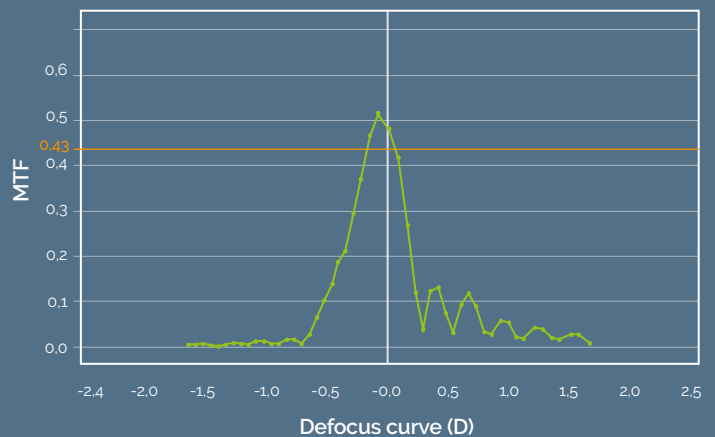

SPECTRAL TRANSMISSION

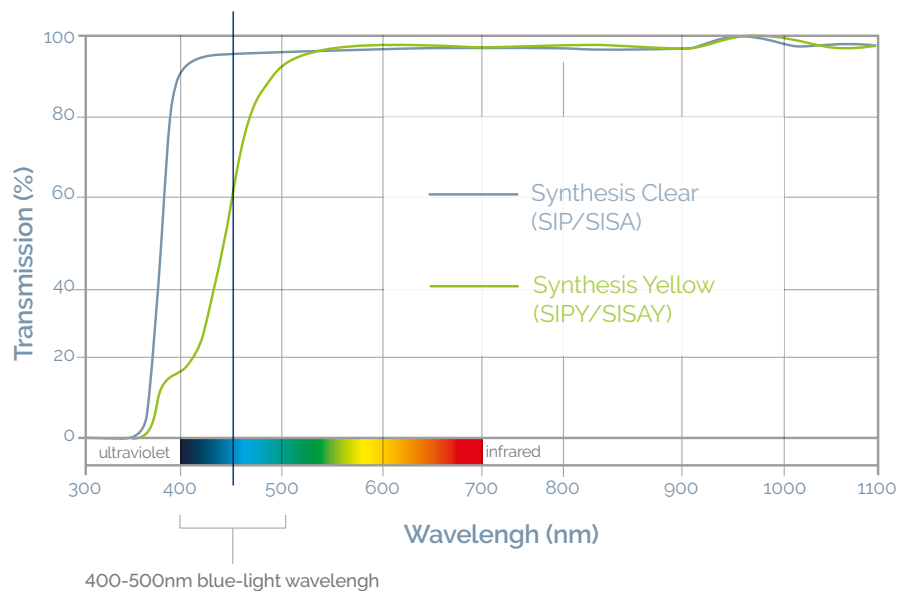

Synthesis Yellow filters 40% of blue light at 450nm wavelength



TECHNICAL SHEET

DESIGNATION	SPECIFICATIONS
Optic diameter	6 mm across the power range (0.0 D to 30.0 D)
Diopter Range	1 D increment between -10.0 D to +10.0 D & +31.0 D to +40.0 D 0.5 D increment between +10.0 D to +30.0 D
Ordering information	Stock items between 0.0 D to 35.0 D On-order (4 weeks lead-time) between -10.0 D to -1.0 D & 36.0 D to 40.0 D
Total diameter	From -10.0 D to 15.0 D: 11.0 mm From 15.5 D to 22.0 D: 10.7 mm From 22.5 D to 30.0 D: 10.5 mm
Design	4 point-fixation, C-Edge+ \leq 5 microns
Angulation	0°, posterior-shift feature
Optic	Aspheric aberration-neutral
Biometry	SRK/T: A = 118.851 Haigis: a0 = 1.2228 (with a1 = 0.4 & a2 = 0.1) Hoffer-Q: pACD = 5.4595 Holladay 1: SF = 1.7064
References	Synthesis <i>CLEAR</i> Preloaded: SIP Synthesis <i>YELLOW</i> Preloaded: SIPY Synthesis <i>CLEAR</i> IOL: SISA Synthesis <i>YELLOW</i> IOL: SISAY

INJECTION SYSTEMS

		Brands	References	Wound-assisted recommended incision	In-the wound recommended incision	In-the-bag recommended incision
	Preloaded	Medicel Accujet™ PRO 1.6	LP604550P (box of 10)	1.6 mm	1.8 mm	2.0 mm
	Cartridge loading	Medicel Viscojet™ BIO 1.8	LP604350C (box of 10)	1.8 mm	2.0 mm	2.2 mm



CUTTING EDGE
sense of innovation

SENSE OF EXCELLENCE

SENSE OF INNOVATION

SENSE OF PARTNERSHIP

SENSE OF ETHICS

CUTTING EDGE HEADQUARTERS

770 rue Alfred Nobel - Immeuble le Nobel
34000 Montpellier - France

CUTTING EDGE MANUFACTURING

580 rue Max Planck - 31670 Labège - France

KESTREL OPHTHALMICS LTD

Unit 9 Cabot Business Village, Holyrood Close, Poole, Dorset, BH17 7BA - United Kingdom
Tel. +44 (0) 1202 696963

CUTTING EDGE SURGICAL

70I2Q Local Avinguda Manuel Girona - 08860 Casteldefels - España
Tel. +34 658 208 691

www.cutting-edge.eu - info@cutting-edge.eu